



PRODUCT MANUAL

A Sterile Processing Users Guide



Safety, Assurance, Innovation, and Quality.

Enabling perioperative professionals to achieve higher standards of performance and enhanced patient safety.



Innovative Sterilization Technologies
7625 Paragon Road, Suite A
Dayton, OH 45459
www.iststerilization.com

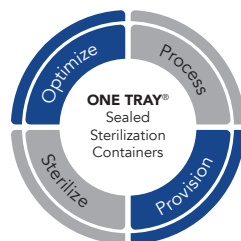


Condoc 593 Rev 20



Proudly Made in the U.S.A.

©2024 All rights reserved. Not to be duplicated without express permission of the author.



CONTENTS

Welcome to ONE TRAY®	3
ONE TRAY® System Recommended Usage & Sterilization Guidelines	4
Shelf Life	4
ONE TRAY® 510k Cleared Performance Testing Summary	5
Preparation	6
Opening Method	7
Disposable® Filter Installation	8
Loading The Container	10
Container Assembly and Security	11
Closing Method	12
Tamper Evident Lock Method	13
Sterilization and Handling	14
Sterile Delivery and Aseptic Presentation	15
Removing Sterile Contents	15
Inspection Checklist	16
Reprocessing Instructions	17
Disassembly	17
Cleaning	18
Manufacturer's Warranty	19
Limits of Use	20
ONE TRAY® Processing Kit Storage Requirements	21

Personnel training and competency is required to perform all phases of processing. Sterilizing equipment, water supply and quality, and practices within a facility all contribute to providing an effective reprocessing system and should be monitored by each facility.

WELCOME TO ONE TRAY®

Thank you for purchasing the ONE TRAY® Sealed Sterilization Container system.

ONE TRAY® Sealed Sterilization Containers are intended to be used to hold temperature tolerant medical devices, surgical supplies, single instruments, multiple instruments or an instrument set for immediate use following flash sterilization (defined as a steam sterilization cycle: Temperature 270°F (132°C), Exposure Time 4 minutes, Cycle Dry Time Not Applicable). ONE TRAY® containers have been validated to maintain the sterility of the contents for 48 hours. This includes sterilization of lumens 3 mm in diameter or larger with lengths of up to 400 mm.

After flash sterilization (defined as a steam sterilization cycle: Temperature 270°F (132°C), Exposure Time 4 minutes, Cycle Dry Time Not Applicable), ONE TRAY® provides for the safe transport and assured delivery for immediate use, and maintenance of sterility of the enclosed devices for 48 hours storage in a sealed container with tamper evident security and load record documentation according to AAMI and AORN guidelines."



For device model numbers, sizes and product accessory information please go to onetray.com/products/onetray for details.

ONE TRAY® SYSTEM RECOMMENDED USAGE & STERILIZATION GUIDELINES

ONE TRAY® is validated to process a twenty five pound (25 lb) gross weight load (single container plus contents) in the following parameters listed in Table 1. This includes sterilization of lumens 3 mm in diameter or larger with lengths of up to 400 mm.

Run loaded sterilizer according to the times and temperatures listed. The following parameters were established during validation testing of ONE TRAY® and performed in an AAMI ST-08 compliant steam sterilizer.

ONE TRAY® is validated for use utilizing the following parameters:

Table 1

Parameter	Dynamic Air Removal/ Pre-Vacuum Steam Sterilization
Exposure Temperature	270°F (132°C)
Exposure Time	4 minutes
Dry Time	Not Applicable
Cool Time	Varies according to load

SHELF LIFE

ONE TRAY® Sealed Sterilization Containers can be used immediately following flash sterilization (defined as a steam sterilization cycle: Temperature 270°F (132°C), Exposure Time 4 minutes, Cycle Dry Time Not Applicable), ONE TRAY® recognizes the hectic nature of the perioperative environment. ONE TRAY® containers have been validated to maintain the sterility of the contents for 48 hours. If devices processed in the ONE TRAY® are not used immediately within the validated 48 hour storage period, they should be reprocessed following normal departmental procedures (i.e. Shelf Life, Storage, etc.)

ONE TRAY® 510K CLEARED PERFORMANCE TESTING SUMMARY

Contents/Configuration	Performance Testing conducted for 510k K052567
Required Accessories	ONE TRAY® container, filters and tamper evident locks
Lumens	Lumens 3mm in diameter or larger with lengths up to 400mm
Shelf Life	ONE TRAY® containers have been validated to maintain the sterility of the contents for 48 hours.
Maximum Weight	Total gross weight, container plus contents not to exceed 25lbs
Device Types	Instruments with hinges, knurled areas and lumens
Device Materials	<u>Metals:</u> Stainless Steel <u>Aluminum Thermoplastics:</u> Polypropylene, Radel <u>Thermosetting Polymers:</u> Silicone



All clinical research and validation testing of ONE TRAY® was performed using the ONE TRAY® container, filters, and tamper evident locks.

Performance testing of the ONE TRAY® was conducted by HIGHPOWER Labs, an ISO 17025 accredited laboratory. All testing was performed without the use of a dry time. The following performance testing / validations were conducted in accordance with the below listed recognized standards:

Performance Testing / Validation:

- Steam Lethality
- Shelf Life
- Sterilant Penetration
- Biocompatibility
- Material Compatibility
- Moisture Sterility
- Package Integrity
- Limit of Reuse

Testing / Validation in accordance to:

- ANSI/ AAMI ST77
- ANSI/ AAMI TIR30
- ANSI/ AAMI ST79
- ANSI/ AAMI ST81
- ANSI/ AAMI ST8
- AAMI 14161
- ANSI/ AAMI TIR12
- ISO 10993
- ISO 11607
- ISO 17665

PREPARATION

1. Thoroughly inspect the container for dents, loose components or damage that may affect performance.
2. Place container on stable and level surface at height that facilitates opening the lid safely.
3. Open the container by removing the lid according to the instructions illustrated in the table on the next page.
4. Confirm the rim of the base is free of burrs and dents.
5. Inspect to ensure all gaskets are free of cracks, tears, imperfections or defects and are completely seated in the applicable channel.
6. Confirm that the lid properly mounts on the base and seals in place.



OPENING METHOD

1	 Lid on base with latch mechanism in locked position	2	 Grasp the latch mechanisms at each end of container
3	 Simultaneously pull outward on the latching mechanisms on each end	4	 Locking lever is disengaged from mounting block on base as the latching mechanisms are raised
5	 Locking lever releases from the mounting block on the base as both latching mechanisms are raised	6	 Latching mechanisms are fully open and lid removed from base



The validated performance of ONE TRAY® is in effect **ONLY** with the exclusive use of the ONE TRAY® container, filters, and tamper evident locks.

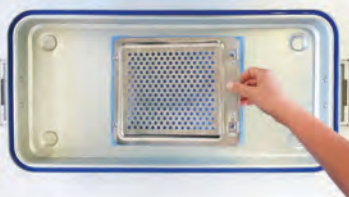
DISPOSABLE FILTER INSTALLATION

1. A single disposable  ONE TRAY® filter is required for each of the perforated areas (one in lid and two in the base).
2. Each disposable  ONE TRAY® filter should be inspected to confirm the absence of imperfections or defects such as puncture holes or tears.
3. Insert the ONE TRAY® filter with the ONE TRAY® name placed face up. The name should be under the filter cover slide bar.
4. Each disposable  ONE TRAY® filter **must** completely cover the grooved channel surrounding the perforated areas.
5. Confirm the notched end of the filter cover is properly oriented with the correct set of retention posts.
6. Insert the notched end of filter cover under retention posts.
7. Position “keyholes” in the filter cover slide-bar over an opposing set of retention posts.
8. Apply downward pressure on the outer edges of the belt loops on the filter cover and push the slide bar completely forward into a fully locked position under the retention posts.
9. Filters and filter covers are interchangeable.



Disposable  filters are for single use only and must be discarded after each processing cycle. Do not use more than one filter in each perforated area that requires a filter.

FILTER INSTALLATION

1	 Use ONE TRAY® Processing Kit	2	 Confirm filter cover seal and grooved channels are undamaged (cracks, tears, burrs, dents, discoloration) and clean
3	 Position the ONE TRAY® filter so that the holes align over the retention posts on the same side as the slide bar on the filter cover. The ONE TRAY® name will align under the slide bar of the filter	4	 Insert notched end of filter cover under appropriate retention post
5	 Position "keyholes" in filter cover slide-bar over the retention posts on the opposite side	6	 Apply downward pressure on the filter cover and push the side-bar into the locked position

LOADING THE CONTAINER

1. Prevent contents from contacting the lid.
2. Avoid obstructing perforated areas and filter covers.
3. A chemical integrator should be used within each assembled content set. Facilities should follow their internal policies and procedures for integrator usage.
4. When loading the container, utilize the ONE TRAY® deck plate or equivalent size containment device for the direct placement of medical devices or instrument organizing trays. Do not set anything directly on the filter covers.



Medical Devices should be prepared and sterilized according to the device manufacturer's written instructions for sterilization cycle exposure parameters (**time and temperature**).

CONTAINER ASSEMBLY & SECURITY

1. Inspect the rim of the base to confirm there are no dents or burrs.
2. Inspect the lid gasket to confirm there are no cracks, tears, imperfections, or defects.
3. Confirm that all filters and filter covers are securely in place.
4. Place the lid on the base evenly and confirm proper fit.
5. Secure the lid by simultaneously pressing down to engage the latching mechanisms at each end.
6. Insert a ONE TRAY® tamper evident lock through the latching mechanisms at both ends of base.

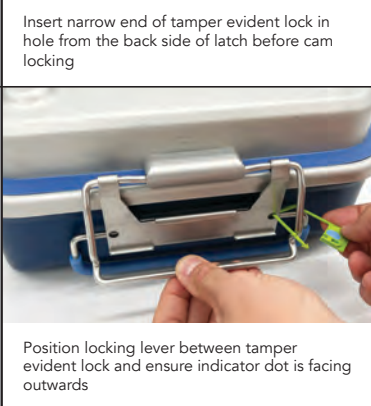
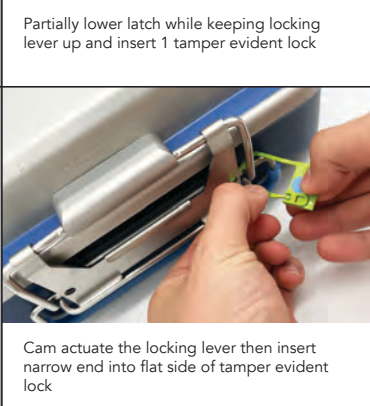


The validated performance of ONE TRAY® is in effect **ONLY** with the exclusive use of the ONE TRAY® container, filters, and tamper evident locks.

CLOSING METHOD

1	 Align the lid with the base	2	 Mount the lid properly onto the base
3	 Engage the latching mechanisms with the mounting block on the base and the locking levers are in the horizontal position	4	 Simultaneously push down on both locking levers
5	 Push down on both locking levers until latches are fully engaged and in the locked position	6	 Container is clamped and sealed when both cover latching mechanisms are in locked position

TAMPER EVIDENT LOCK METHOD

1		2	
3		4	
5		6	

Insert narrow end of tamper evident lock in hole from the back side of latch before cam locking

Partially lower latch while keeping locking lever up and insert 1 tamper evident lock

Position locking lever between tamper evident lock and ensure indicator dot is facing outwards

Cam actuate the locking lever then insert narrow end into flat side of tamper evident lock

Push narrow end through opening until tamper evident locks

Tamper evident lock placed on latch at both ends of container

STERILIZATION & HANDLING

1. Place a single ONE TRAY® container in a level horizontal position on the shelf of the sterilizer or sterilizer rack.
2. Confirm that all perforated areas are not obstructed.



DO NOT STACK ONE TRAY® CONTAINERS DURING STERILIZATION

3. ONE TRAY® is validated for processing a total gross weight (container plus contents) not to exceed 25 lbs. This includes sterilization of lumens of 3 mm in diameter or larger with lengths of up to 400 mm.
ONE TRAY® is validated for use in the following sterilization cycle:

Pre-vacuum

Temperature 270°F (132°C), Exposure Time 4 minutes, Cycle Dry Time Not Applicable



Following sterilization, avoid unprotected contact with the hot container by using insulated handles and/or gloves to lift and transport.

4. Retained moisture is normal and expected and will vary by load composition. This moisture does not compromise the sterility of contents during storage.



Items sterilized in a steam sterilization cycle: Temperature 270°F (132°C), Exposure Time 4 minutes, Cycle Dry Time Not Applicable) are validated for a 48 hour storage period.

STERILE DELIVERY & ASEPTIC PRESENTATION

1. Place container on stable and level surface at height that facilitates aseptic opening.
2. Check to ensure tamper evident locks are intact by gently tugging on each lock.
3. Check to ensure the external indicator dot has turned the appropriate color. Upon exposure to high temperature steam, the indicator will transition from blue to dark. The transition color may vary depending on the load configuration, length and condition of exposure. A color transition from blue to a shade of black/grey provides indication of exposure to high temperature steam.
4. Before opening, inspect the container and confirm the presence of all three filters.
5. Break the tamper evident locks by twisting the lock until the lock breaks. Remove the lock.



Failing to follow ONE TRAY® lock breaking procedures could result in lock pieces transferring onto the sterile field.

6. Place fingers under lid locking lever at both ends of container.
7. Aseptically lift lid off the base of the container and inspect for proper placement of all filters and installation of filter covers. Ensure filter covers are securely in place.

REMOVING STERILE CONTENTS

1. Before removing the contents within the container, trained personnel should verify internal sterilization integrator per facility policies and procedures.
2. Securely grasp the contents, or modular organizing tray handles and remove in an aseptic manner without touching the outside of the container.
3. Before placing contents on the sterile field, each disposable ONE TRAY® filter should be inspected to confirm the absence of imperfections or defects such as puncture holes or tears.

INSPECTION CHECK LIST



Contents should be considered **non-sterile** if any of the following conditions are present:

- A tamper evident lock is not intact or is missing from either latching mechanism.
- A filter cover is dislodged and not securely attached.
- A filter does not completely cover the channel around the vented area.
- Integrity of the filters has been compromised. Each disposable ONE TRAY® filter should be inspected to confirm the absence of imperfections or defects such as puncture holes or tears.
- Gasket is either damaged or separated from the lid channel.
- Sterilization cycle is interrupted or aborted (doesn't complete required steam sterilization parameter).

If you have questions about the proper functionality of the ONE TRAY® Sealed Container System, contact your ONE TRAY® representative.



The validated performance of ONE TRAY® is in effect **ONLY** with the exclusive use of the ONE TRAY® container, filters, and tamper evident locks.

REPROCESSING INSTRUCTIONS

Follow facility's policies, procedures, and AAMI ST79 recommended guidelines for the transportation of soiled instruments and containers.

It is recommended that containers be reprocessed as soon as is reasonably practical following use. Containers should be reprocessed in the completely open and disassembled configuration.

Container, lids and baskets that may not be used or needed right away should be decontaminated and cleaned prior to storage.

The ONE TRAY® Container System should be stored neatly, either assembled or unassembled, in a dry, clean area.

DISASSEMBLY

1. Remove and discard all disposable components (filters, tamper evident locks).
2. Inspect the container lid and base for dents or damage that may affect performance. Confirm the rim of the container base is free of burrs and dents.
3. Inspect to ensure the lid gaskets are free of cracks, tears and is properly seated in the lid channels.
4. Inspect filter cover seal and grooved channels. Ensure these features are undamaged (cracks, tears, burrs, dents, discoloration) and clean.
5. Confirm the lid properly interfaces and engages on container base.

CLEANING

1. ONE TRAY® containers and components may be safely reprocessed in a mechanical washer/decontaminator, cart washer, or processed manually.
2. Mechanical or automatic reprocessing should follow the directions of the equipment manufacturer. The ONE TRAY® lid and base should be positioned to prevent water collection.
3. Use reprocessing agents within a pH range of 6.5 to 8.5.
 - 3.1 If white residue is observed within the container, this may have been caused by the use of reprocessing agents outside the recommended pH range.



All reprocessing agents **must** be thoroughly rinsed off. AAMI TIR 34 compliant water is recommended for the final rinse to avoid discoloration or damage resulting from minerals found in utility water.



Use of reprocessing agents below or above the recommended pH range could permanently damage the protective finish of the container and will result in warranty exclusion.

4. **NEVER** clean the container with abrasives or wire brushes.



Use of abrasive materials will permanently damage the protective anodized finish of the container. Use of abrasive materials will result in warranty exclusion.

5. All components should be secured or enclosed inside a rack or basket with no protrusions to prevent damage during reprocessing.
6. Confirm that all components are cool before handling.
7. Ensure container and contents are dry prior to sterilization. If necessary, dry container and components with a clean lint-free towel.

MANUFACTURER'S WARRANTY

ONE TRAY® Sealed Sterilization Containers are guaranteed for the life of the product to be free of functional defects in workmanship and materials when used normally as recommended for their intended purpose in conjunction with ONE TRAY® filters, tamper evident locks, and deck plate.

ONE TRAY® containers must be inspected per the inspection checklist prior to each usage for signs of misuse or mishandling (i.e. dents, missing parts, damaged pieces, misalignment etc.) and sent in for warranty claims.

All products must be thoroughly decontaminated and cleaned before being returned for warranty claims.

Any product determined to be defective after normal usage will be repaired or replaced with no charge to the customer at the sole discretion of Innovative Sterilization Technologies. A restocking fee will be charged on returns.

Used product may not be returned for credit (used is considered any box that has been cut opened). Any and all non-used returns **must** receive prior authorization from Innovative Sterilization Technologies.

This warranty is valid only to the original purchaser.

The following exclusions apply to this warranty that include, but are not limited to:

- Conditions resulting from: negligence, misuse and improper reprocessing, handling, closing or opening.
- Damage from excessive force or pressure.
- Not using reprocessing agents within a pH range of 6.5 to 8.5
- Modification(s) to the container.
- Effects from fire, flood or other unpredictable event not under the control of Innovative Sterilization Technologies.



ONLY IST trained technicians are authorized to repair ONE TRAY® Sterilization Container System. Using a non-IST repair technician to repair containers will void the ONE TRAY® warranty on the container and will void any of the validation/performance testing associated with ONE TRAY® containers.

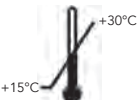
LIMITS OF USE

Repeated reprocessing has minimal effect on the ONE TRAY® device. The useful life depends on many factors including the method and duration of each use, and the handling between uses. Careful inspection before use is the best method of determining the end of serviceable life.

Do not use devices that show evidence of damage and wear. Evidence of damage and wear may include but is not limited to corrosion (i.e. rust, pitting), discoloration, excessive scratches, flaking, wear and cracks.

Ref #	Symbol	Title of Symbol	Description of Symbol
5.1.1		Manufacturer	Indicates the medical device manufacturer.
5.4.2		Do not re-use	Indicates a medical device that is intended for one use, or for use on a single patient during a single procedure.
5.4.3		Consult Instructions for Use	Indicates the need for the user to consult the instructions for use.
5.4.4		Caution	Indicates the need for the user to consult the instructions for use for important cautionary information such as warnings and precautions that cannot, for a variety of reasons, be presented on the medical device itself.

ONE TRAY® PROCESSING KIT STORAGE REQUIREMENTS

	15°C to 30°C		Keep away from direct exposure to sunlight and fluorescent light
	20% to 70% Relative Humidity		Keep Dry
	Reference UDI device label attached to the ONE TRAY® processing kit box for expiration date.		
	Keep away from hydrogen peroxide and other sterilants. Storage on top of or near a heat source should be avoided. Do not use after expiration date.		

- Industry Definitions can be found on the ONE TRAY® Website at: <https://onetray.com/definitions>
- Vent to volume ratios for the ONE TRAY® container can be found on the ONE TRAY® website at: <https://onetray.com/onetrayspecs>

CONTACT

If you have questions regarding the **ONE TRAY®** warranty, contact

info@onetray.com or
call 937.619.0138.