

ONE TRAY®
Sealed Sterilization Container

STERILIZATION MADE **SIMPLE.**



ONE TRAY®
Sealed Sterilization Container

ONETRAY.COM

ONE TRAY®

Sealed Sterilization Container



- Faster Instrument Turnaround
- Increased OR Throughput
- Validated Sterility Maintenance

ONE TRAY® is an American made, lifetime warranted sealed sterilization container offering efficient instrument reprocessing in a fraction of the time it takes sterile wrap or other rigid containers.

REPROCESSING TIME SAVINGS EXAMPLE

ONE TRAY®

Sealed Sterilization Container



Ready for the OR **2 hrs+** sooner!

Watch this LIVE Knee Case
Timelapse video showing
ONE TRAY® vs Sterile Wrap.



OTHER PACKAGING METHODS




APPROX. 2.5+ HOURS

ADDITIONAL DRY AND COOL TIME REQUIRED

ONE TRAY® FDA CLEARANCES


IFU K052567

DRY TIME NOT REQUIRED, 48 HR STORAGE


Exposure	270°F (132°C)	
Exposure	4 minutes	
Dry Time	Not Required	
Cool Time	Not Required	
Validated Storage Period	48 hrs	

IFU K250029

15 MIN MINIMUM DRY TIME, 365 DAY STORAGE

Exposure	270°F (132°C)	
Exposure	4 minutes	
Dry Time	15 min minimum	
Cool Time	Varies according to load	
Event Related Shelf Life/ Storage Period	365 days	

PROCESSING KITS



The **ONE TRAY®** processing kits include single use disposable  filters that permits the penetration of steam during the sterilization cycle and serves as a bacterial and fluid barrier at the conclusion of the sterilization cycle.

Kits include*:

- 3 Filters - disposable hydrophobic SMS Filters
- 1 Chemical Integrator
- 4 Locks with indicator dots - (2) Green for use with IFU K052567 and (2) Orange for use with IFU K250029
- 2 Load Cards - for easy identification

*Locks may come packaged separately

The validated performance of **ONE TRAY®** is in effect **ONLY** with the exclusive use of **ONE TRAY®** components, filters, and tamper evident locks.

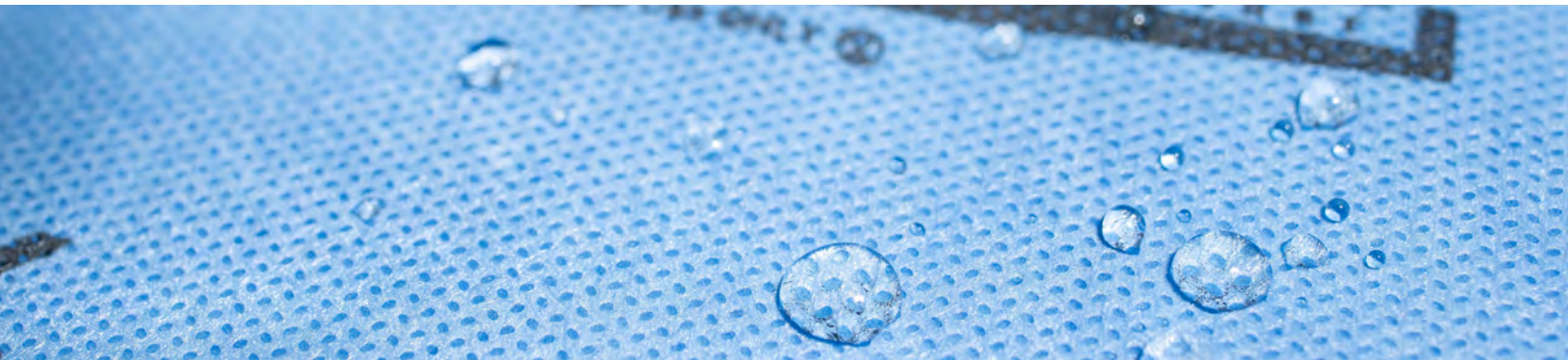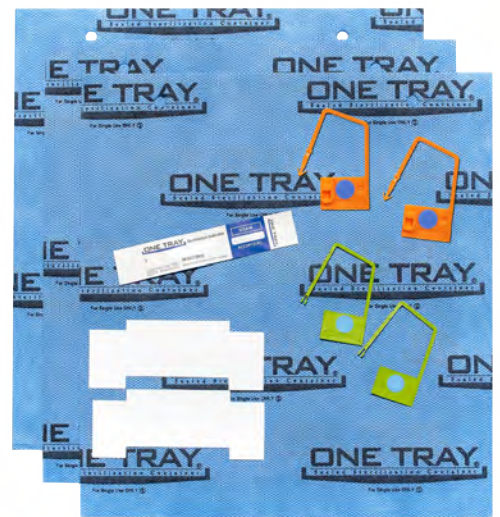

TRAY MODELS AND DESCRIPTIONS

Container Catalog #	M2104 M2104DT	M2106 M2106DT	M2108 M2108DT	M2404 M2404DT	M2406 M2406DT	M2408 M2408DT
K052567/ K250029	21"l x 12"w x 4"h with 21" deck plate	21"l x 12"w x 6"h with 21" deck plate	21"l x 12"w x 8"h with 21" deck plate	24"l x 12"w x 4"h with 24" deck plate	24"l x 12"w x 6"h with 24" deck plate	24"l x 12"w x 8"h with 24" deck plate



M2104

M2404

M2106

M2406

M2108

M2408



Filter Covers use a patented **SterObex™** barrier providing sealing capabilities for sterility maintenance

Our mission is to redefine sterilization through efficient, scalable solutions that uphold the highest standards of patient care.



EVERY MINUTE MATTERS. LET'S MAXIMIZE YOURS.



Scan here and let's talk about optimizing your surgical workflow.

937.619.0138
info@onetray.com



Innovative Sterilization Technologies
7625 Paragon Road, Suite A, Dayton, OH 45459
937.619.0138 | www.onetray.com



Reference the **ONE TRAY®** IFU's for usage instructions.



Made in the USA
Condoc 336 Rev03